


Key Features

- Large platform rockers – ideal for multiple users
- Gentle rocking action
- Choice of two models:
 - 3D gyratory action SSL3
 - See-saw wave action SSL4
- Digital speed control and built-in timer
- Optional tier system available to increase capacity



SSL3

SSL3 and SSL4

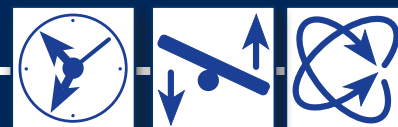

Rockers, lab scale,

These rockers have large platforms able to accommodate a number of samples, ideal for a busy lab. They are very quiet in operation and designed to be on continuously. Two models available:

Model **SSL3** provides a 3D gyratory motion, ideal for low foaming agitation, DNA extractions, staining and de-staining procedures etc. The angle of tilt can be moved to any position by hand to optimise mixing of vessels.

Model **SSL4** has a see-saw rocking action that creates a wave motion within vessels such as culture flasks, petri dishes etc.

Rockers are often used in conjunction with incubators and environmental chambers. Both units can be used in temperatures up to 40°C and humidity up to 80%. Both models are supplied with a non-slip mat and have digital selection of both speed and time making them very easy to operate. An accessory tier system can be fitted in seconds, without the need for tools, tripling the available space for samples. Each tier is securely held in place by magnets.

Technical Specification

| | SSL3 | SSL4 |
|---------------------------------|-----------------|-----------------|
| Rocking action | 3D gyratory | See-saw wave |
| Platform dimensions, mm (w x l) | 355 x 355 | 355 x 355 |
| Speed range | 5 to 70rpm | 5 to 70rpm |
| Angle of tilt | 3 to 12° | 7° |
| Angle adjustable | Yes by hand | No |
| Maximum load, kg | 10 | 10 |
| Dimensions, mm (w x d x h) | 360 x 420 x 170 | 360 x 420 x 160 |
| Tier height, mm | 125/tier | 125/tier |
| Operational temp. range | +4 to +40°C | +4 to +40°C |
| Maximum permissible humidity | 80% | 80% |
| Net weight, kg | 10 | 10 |
| Electrical supply | 230V, 50Hz, 50W | 230V, 50Hz, 50W |
| IP Rating | 30 | 30 |

Ordering Information

| Model | Description |
|--------|---|
| SSL3 | Rocker, gyratory, lab scale |
| SSL4 | Rocker, see-saw, lab scale |
| SSL3/1 | Tier system (2 platforms plus 8 bars) (fits both models) |



SSL4